

Surface Activity Monitor


Model: 7001

Application

The Surface Activity Monitor quickly, conveniently and directly measures fixed and removable tritium activity on flat surfaces and can measure removable tritium on all surfaces. This instrument is portable, can have its calibration easily checked, produces no waste, and is much quicker and more convenient than scintillation testing for tritium on surfaces.

Features

- Measures removable tritium on smear paper.
- Measures total tritium on flat conducting surfaces.
- No P10 gas required.
- Simple and easy to use.
- No waste.
- Portable.
- Optional calibrated source is available.

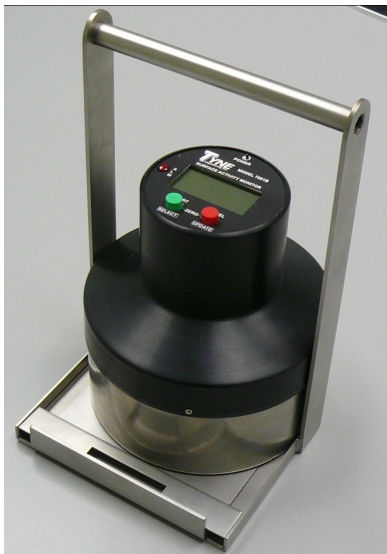

Description

Smearing or wiping, is the most common and reliable method of determining tritium surface contamination in pipes, on the surfaces of drums, glove boxes, laboratory benches, or any other devices that may have tritium contamination on surfaces. A smear is normally measured by dissolving the smear paper in a "cocktail" putting it into a vial and measuring disintegrations

in a scintillation counter. There are subsequent problems associated with storage of the active liquid and cleaning the apparatus. The scintillation counter is known to be an accurate instrument, but since making a smear sample is less accurate, it is the efficiency of creating the smear that controls the overall accuracy of the process.

Tyne has developed an instrument that overall is as accurate as the scintillation counter for reading smears. With the SAM monitor, it takes about 2 minutes to arrive at a smear activity measurement. In comparison, the time for a scintillation measurement may be hours or days depending on the laboratory facilities available.




Read-out and controls

The SAM, which is battery operated, may be placed close to an operator on a laboratory bench for maximum convenience when taking smears. The procedure is simply to open a drawer in the SAM, place in the smear paper, wait 2 minutes and read the measurement.

Measurements are displayed on a large LCD screen on the top of the instrument. Various configuration settings such as alarm settings, measurement units, etc. can be configured using two push buttons on the top panel.

Optionally Tyne can provide an amorphous tritium source of known activity which can be used daily to confirm the calibration of the instrument. If a particular reading is important, a calibration check can easily be made both before and after the reading. The calibration-checking source is mounted in a plate which slides into the same drawer used for the smear paper.

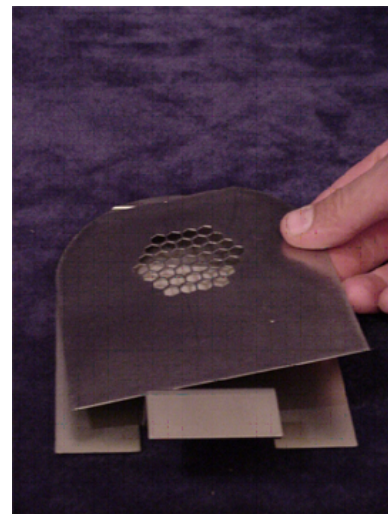

mounting plate for smear paper

The SAM can also be used directly, without smear paper on flat conducting surfaces such as a stainless steel drum or glove box, or a laboratory bench to check for tritium on surfaces. A reading taken directly, measures not only the removable activity which can be picked up in a smear paper, but also the non-removable near surface contamination which would otherwise not be recorded by the smear.

Specifications

Model	Low Range (Tyne P/N - 00047602)
Range	0-3552 kDPM/dm ² (10-6000 nCi/dm ²)
Detection Sensitivity	+/-16.65 kDPM/dm ² (+/- 7.5 nCi/dm ²)
Response Time	<2 minutes
Measurement Area	100 cm ²
Display	128x64 Graphic LCD Display
Accuracy	+/- 4.44 kDPM/dm ² (0-100 kDPM/dm ²), +/-10% above 100kDPM/dm ²
Process/Ambient Temperature	15°C to 50°C
Relative Humidity	5 to 65%
Air Pressure	80-120kPa
Tritium Wetted Parts	316L Stainless Steel, Teflon
Dimensions	150 mm dia x 178 mm high x 254 mm wide
Weight	2.7 kg
QA	ISO 9001:2000
Calibration	An optional calibrated source comprising amorphous silicon tritium emitter can be provided. Please call for availability and pricing
Source strength	22- 2,200 KDPM (10-1000nCi/dm ²)
Units	kDPM/dm ² , nCi/dm ² , or Bq/cm ²
Alarm level setting	Full range.